





Rec'd PCT/PTO 17 MAY 2005  
GB 03/4965



INVESTOR IN PEOPLE

The Patent Office  
Concept House  
Cardiff Road  
Newport  
South Wales  
NP10 8QQ

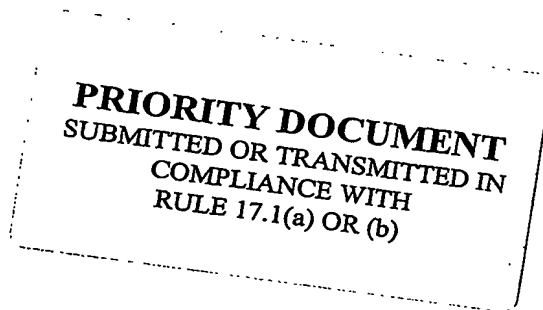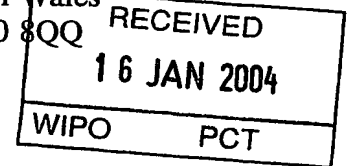


I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

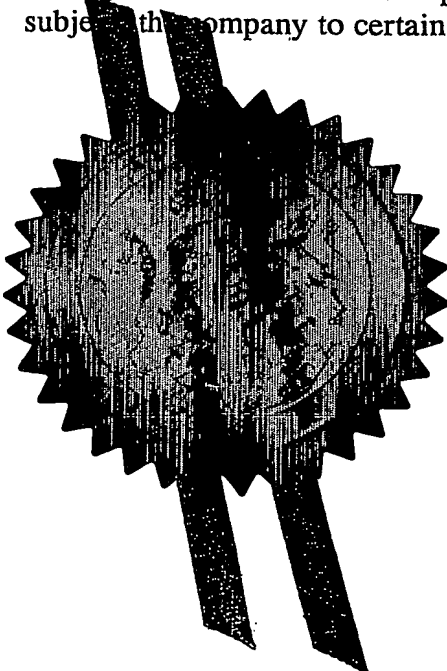

Signed

*Stephen Hordley*

Dated 19 December 2003

**BEST AVAILABLE COPY**

20 NOV 2002

NEWPORT



20NOV02 09:48:61-1 C58534  
P01/7700 0.00-0227043.7

The Patent Office

Cardiff Road  
Newport  
South Wales  
NP10 8QQ

# Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

1. Your reference

JA/bk/2591GB

2. Patent application number

(The Patent Office will fill in this part)

0227043.7

3. Full name, address and postcode of the or of each applicant (underline all surnames)

Smith & Nephew, plc  
15 Adam Street  
LONDON WC2N 6LA

Patents ADP number (if you know it)

03969284006 ✓

If the applicant is a corporate body, give the country/state of its incorporation

4. Title of the invention

Angiogenic Medical Cyanoacrylate  
Adhesive.

5. Name of your agent (if you have one)

J. D. Hobbs

"Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)

Smith & Nephew Group Research Centre  
York Science Park, Heslington,  
York YO10 5DF  
United Kingdom

Patents ADP number (if you know it)

1116615003

6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (if you know it) the or each application number

Country

Priority application number  
(if you know it)

Date of filing  
(day / month / year)

7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application

Number of earlier application

Date of filing  
(day / month / year)

8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (Answer 'Yes' if:

YES

- a) any applicant named in part 3 is not an inventor, or
  - b) there is an inventor who is not named as an applicant, or
  - c) any named applicant is a corporate body.
- See note (d))

## Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form. Do not count copies of the same document

Continuation sheets of this form

Description

9

Claim(s)

Abstract

Drawing(s)

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (Patents Form 7/77)

Request for preliminary examination and search (Patents Form 9/77)

Request for substantive examination (Patents Form 10/77)

Any other documents (please specify)

11.

I/We request the grant of a patent on the basis of this application.

Signature

Date 19/11/02

JOHN HOBBS

12. Name and daytime telephone number of person to contact in the United Kingdom

John Hobbs 01904 824050

### Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

### Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 08459 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

## ANGIOGENIC MEDICAL CYANOACRYLATE ADHESIVE

The invention relates to adhesives and angiogenesis, in particular the use of novel angiogenic adhesives in surgery.

5

Adhesion of tissue is an integral part of all surgical procedures, including closure of skin wounds, reconstruction of nerve ruptures, re-attachment of transplanted tissue, sealing of blood vessels, treatment of pneumothorax and fistulas, support of vascular and  
10 intestinal anastomoses, treatment of chondral- and osteochondral defects, fracture healing, treatment of meniscal tears and ruptured ligaments, repair of tendon damage or muscle damage and attachment of implanted biomaterials and tissue engineered devices.

15

The fundamental aim of all tissue adhesives is to hold tissue together for long enough to allow a natural biological repair. Biological repair typically involves the activation of cells in the tissue to a repair mode and, critically, the stimulation of angiogenesis to provide repair cells, nutrition and oxygen to the activated cells.

20

There is extensive literature describing the use of bioactives to stimulate tissue repair. However, stimulation of tissue repair mechanisms alone, such as angiogenesis, using biological actives will not be enough to heal a large defect between tissues. It is  
25 essential that the two tissues to be integrated are held in close apposition at the macroscopic level in order to allow the biological mechanisms to bond the tissue at the microscopic level.

To date tissue adhesion is mainly done using mechanical  
30 fastening techniques such as suturing or stapling. There are occasions, however, when the application of a biological glue or adhesive would be beneficial. For example, sutures are inappropriate for cartilage repair as these cause non-healing defects to form in the cartilage where they are placed. Suturing of internal  
35 tissue and organs is also slow and technically difficult compared to application of an adhesive. Other tissues that may require adhesion such as bone or certain implants may be too hard for sutures or

staples whilst other soft tissues may be too fragile for the suture or staple to hold under tension.

5       Consequently, adhesives have been developed for biological applications, including biological adhesives such as fibrin and synthetic adhesives such as cyanoacrylates.

10       Biological adhesives that utilise naturally occurring adhesive processes such as the blood coagulation cascade (fibrin) have a number of advantages. They are readily accepted by the body and break down completely to allow a full biological repair. However, the bonding strength of such adhesives is well below the levels required for many applications, including all those where the bonded tissue is under any significant tension.

15       A number of synthetic adhesives have been manufactured for industrial and consumer use. Some of these, including cyanoacrylates, have been used to glue biological tissues.

20       The advantage of using cyanoacrylates is that they form an extremely strong bond between tissues. However, they have not replaced the use of other fixation devices because the cyanoacrylate acts as a barrier to biological repair.

25       It would therefore be desirable to produce a cyanoacrylate adhesive that actively stimulated tissue repair.

30       It is a well established fact that, in most tissues, stimulation of angiogenesis results in the acceleration of tissue repair whilst inhibition of angiogenesis can result in ischaemia and tissue death.

35       There is published literature on the delivery of angiogenic factors to stimulate tissue repair. WO 97/16176, WO 01/03607 and US 6,152,141 describe the release of angiogenic factors to accelerate blood vessel repair, EP 0,295,721 describes the promotion of meniscal healing with angiogenic factors, whilst EP 0,530,804 describes the use of angiogenic materials to promote

cartilage and bone healing. However, none of these prior art teach the incorporation of angiogenic factors into cyanoacrylate adhesives.

5 Butyric acid is a potent angiogenic agent and has been used as an angiogenic factor for the treatment of burns, wounds and bone fractures. Butyric acid, also known as butanoic acid, is a four carbon fatty acid. The literature suggest that the local release of 10-1000ng of butyric acid is sufficient to achieve the desired angiogenic effect. However, butyric acid is known to be removed rapidly from the body  
10 and therefore for therapeutic angiogenic applications it has been suggested that it be incorporated into a sustained release delivery vehicle. A lipid angiogenic factor has been isolated from omentum (Catsimpoolas et al., 1984, JAMA 252:2034-2036). The angiogenic factor was found to be monobutyryl (Wilkinson et al., 1991, J. Biol.  
15 Chem. 266:16886-16891).

Monobutyryl can be considered to be a prodrug of butyric acid. Other prodrugs include tributyrin. Tributyrin can be hydrolysed to release butyric acid (Chen *et al*, 1994, Cancer Research 54, 3494 –  
20 3499, Bohmig *et al*, 1999, Transplant Immunology, 7, 221-227). Tributyrin has been proposed for use in anti-cancer therapies where it is desirable to inhibit angiogenesis, it has not been considered as an angiogenic drug.

25 There is prior art describing the use of cyanoacrylate adhesives in tissue repair applications (Barley et al., U.S. Pat. No: 6,342,213, methods for treating non-suturable wounds by use of cyanoacrylate adhesives; Hyon et al., U.S. Pat. N<sup>o</sup>: 6,316,523, adhesive composition for surgical use; Shalaby, U.S. Pat. N<sup>o</sup>: 6,299,631,  
30 polyester/cyanoacrylate tissue adhesive formulations; Kotzev, U.S. Pat. N<sup>o</sup>: 6,224,622, bioabsorbable cyanoacrylate tissue adhesives). However, none of these patents disclose the addition of an angiogenic component to the cyanoacrylate.

35 There is literature describing the incorporation of active molecules into tissue sealants (MacPhee et al., U.S. Pat. N<sup>o</sup> 6,117,425, Supplemented and unsupplemented tissue sealants,

method of their production and use), but this patent makes no mention of cyanoacrylates, despite providing long lists of other frequently used tissue sealants.

5           Simple active molecules such as Iodine have been incorporated into cyanoacrylates (Askill et al., U.S. Pat. N°: 6,214,332, methods for closing suturable wounds by use of cyanoacrylate ester compositions comprising an antimicrobial agent), but angiogenic agents are usually proteins or chemically  
10   active nucleophiles that will cure the cyanoacrylate prematurely, rendering it useless as a tissue adhesive. Thus there is prejudice in the art that suggests that it would not be possible to incorporate an angiogenic agent into a cyanoacrylate.

15           It is an objective of the present invention to provide a cyanoacrylate adhesive for use in biological applications that releases a biologically active angiogenic agent.

20           We have discovered, surprisingly, that tributyrin and some related molecules (butyric acid prodrugs) are also capable of stimulating angiogenesis and do not cause cyanoacrylate to prematurely cure.

25           Accordingly, to the present invention there is provided a tissue adhesive comprising a cyanoacrylate in combination with an angiogenic factor, which is releasable in amounts that will cause a pharmacological effect.

30           The cyanoacrylate adhesive will typically be selected from the group consisting of alkyl 2-cyanoacrylate, alkenyl 2-cyanoacrylate, alkoxyalkyl 2-cyanoacrylate, and carbalkoxyalkyl 2-cyanoacrylate, wherein the alkyl group of said one or more cyanoacrylates has 1 to 16 carbon atoms.

35           The cyanoacrylate will preferably be selected from the group consisting of methyl 2-cyanoacrylate, ethyl 2-cyanoacrylate, n-propyl 2-cyanoacrylate, iso-propyl 2-cyanoacrylate, n-butyl 2-cyanoacrylate,

iso-butyl 2-cyanoacrylate, hexyl 2-cyanoacrylate, n-octyl 2-cyanoacrylate, 2-octyl 2-cyanoacrylate, 2-methoxyethyl 2-cyanoacrylate, 2-ethoxyethyl 2-cyanoacrylate and 2-propoxyethyl 2-cyanoacrylate.

5

In a first embodiment of this invention the angiogenic factor is butyric acid or a derivative or precursor thereof.

10 The angiogenic factors may include:

- Butyric acid (butanoic acid,  $C_4H_8O_2$ ) and butyric acid salts, including sodium, potassium, calcium, ammonium and lithium salts
- 15 • Butyric acid derivatives and polymers containing butyric acid residues
- $\alpha$ -monobutylin (1-glycerol butyrate; 1-(2,3 dihydroxypropyl)butanoate;  $C_7H_{14}O_4$ )
- $\alpha$ -dibutylin (1,3-glyceroldibutyrate; 1,3-(2 hydroxypropyl)dibutanoate;  $C_{11}H_{20}O_5$ )
- 20 •  $\beta$ -dibutylin (1,3-glyceroldibutyrate; 1,2-(3 hydroxypropyl)dibutanoate;  $C_{11}H_{20}O_5$ )
- tributyrin (glycerol tributyrate; 1,2,3-(propyl)tributanoate;  $C_{15}H_{26}O_6$ )
- 25 • hydroxybutyrate and polymers containing hydroxybutyric acid residues

The angiogenic factor is added to the cyanoacrylate in such proportions as to result in an adhesive strength, aptly of not less  
 30 than 0.05Mpa, preferably at least 0.2MPa and more preferably at least 0.5 Mpa. Typically the resultant adhesive strength should range from 0.05 to 0.8 Mpa.

Once cured, the cyanoacrylate will aptly release at least 1ng/ml  
 35 of the angiogenic factor. Suitably the cyanoacrylate will release less than 100 $\mu$ g of angiogenic factor, although preferably it will



release less than 10 $\mu$ g and more preferably less than 1 $\mu$ g of the angiogenic factor.

5 The invention will be illustrated by references to the following examples and the accompanying drawings in which Figure 1 shows that 5% to 50% (w/w) tributyrin can be added to cyanoacrylate without an unacceptable loss of adhesive property. Figure 2 shows that tributyrin is released from a 5% w/w tributyrin cyanoacrylate sample over a 7 day period.

10

Example 1.

0.5g of the sterile tributyrin was added to cyanoacrylate (9.5g) to produce a 5% (w/w) blend in a sterile plastic universal tube. This  
15 was mixed for 2 minutes to ensure a homogenous blend. 40 $\mu$ l aliquots of the blended adhesive were applied to freshly cut surfaces of costal cartilage. The two surfaces were held together for 1 minute allowing fixation and curing was allowed to continue for an additional hour in aqueous conditions. The bond strength was tested using a  
20 Nene MC3000 machine, where the applied force required to separate the bonded pieces of cartilage was recorded. The data showed that 5% to 50% (w/w) tributyrin can be added to cyanoacrylate without an unacceptable loss of adhesive property.

25 Example 2.

0.5g of the sterile tributyrin was added to cyanoacrylate (9.5g) to produce a 5% (w/w) blend in a sterile plastic universal tube. This was mixed for 2 minutes to ensure a homogenous blend. Small  
30 disks of the blended adhesive were cast onto a basic solution (dilute triethylamine 0.1% aq). The disks were removed, washed briefly and dried. The disks were then placed into 2 ml water and stored at 50°C for 1, 5 or 7 days with continuous agitation. The water was removed then added to 2 mls dichloromethane in order to extract any released tributyrin. The dichloromethane layer was analysed  
35 using gas chromatography and the amount of tributyrin recorded. The disks were dried and then placed into fresh water at 50°C for an additional day. The amount of tributyrin released into the fresh water

was measured as described for the first samples. The data showed that tributyrin is released from a 5% w/w tributyrin cyanoacrylate sample over a 7 day period.

Figure 1. Effect of addition of tributyrin to cyanoacrylate adhesive strength


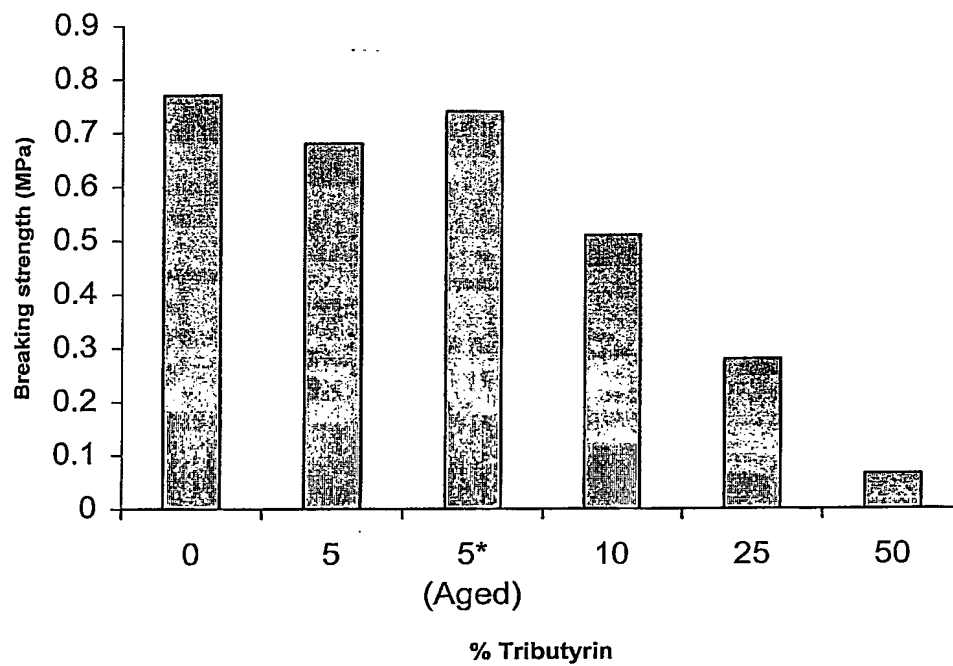


Figure 2. Tributyrin release from 5% tributyrin cyanoacrylate kept at 50°C in water

5

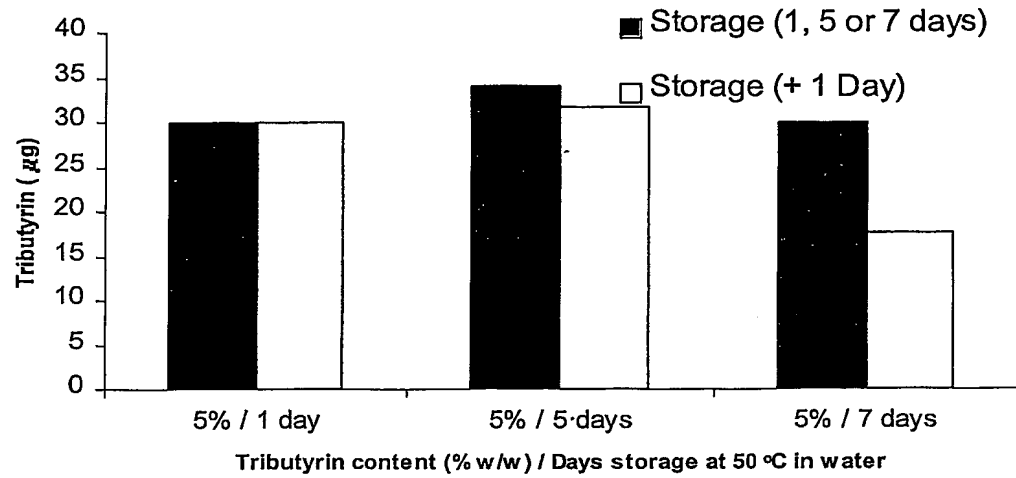

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**